















## Specifications

For the Kiddo Classic the following data are applicable.

	Seat width 22 - 43 cm		Overall length from 72 cm
	Seat depth 20 - 45 cm		Overall width from 52 cm
	Seat height 40 - 58 cm		Back angle 10 degrees
	Leg length 16 - 45 cm		Camber 7 degrees
	Back height 21 - 48 cm		Maximum load 75 kg.
	Overall weight from 17,4 kg included complete seat cushion and back cushion		Rear wheel size 22" 24" 26"
	Weight without rear wheels from 12,4 kg included complete seat cushion and back cushion		Front wheel size 5" 6"

## Standardization

Maximum user weight is 75kg. The Kiddo Classic has been tested and approved by TÜV Rheinland for NEN-EN 12182 and NEN-EN 12183 and crash test passed for occupied transportation: ISO7176 -19. Veldink4kids is affiliated with Nationaal Keurmerk Hulpmiddelen (NL)



# MURRAY'S

## Classic


Kiddo by Veldink

The Kiddo Classic is a growth adjustable wheelchair with strong, rigid frame for everyday use. Customized to the perception of the young user.

Do you want to contact us?

 Block C Airtown Park,  
Airtown Rd, Tallaght,  
Dublin 24


 Eire Tel: 01 855 5696  
UK Tel: 028 91 877951





 info@murrays.ie  
sales@murraysmedical.co.uk

Do you want to follow us?

 [www.facebook.com/murrays.ie](http://www.facebook.com/murrays.ie)  
[www.facebook.com/murraysuk](http://www.facebook.com/murraysuk)

 [www.twitter.com/murraysmed](http://www.twitter.com/murraysmed)  
[www.twitter.com/murraysmeduk](http://www.twitter.com/murraysmeduk)

 [www.murrays.ie](http://www.murrays.ie)  
[www.murraysmedical.co.uk](http://www.murraysmedical.co.uk)

-  Robust, light moving and manoeuvrable wheelchair
-  Growth adjustable in 3 dimensions
-  Wide range of extra options and personal products
-  Perfectly adjustable to each user's therapeutic and physical demands

 Veldink4kids  
no limits

# Kiddo Classic

The child was the central point of interest for developing the Kiddo Classic. This 4-wheel children's wheelchair is custom made and is a combination of design and functionality.

No child is the same and therefore no Classic is the same. Each Kiddo Classic is built on the needs of the user. For the Classic is a wide variation in measurements possible. The wheelchair can grow with the user. The chair is easy to adjust. See size chart for details.

*Excellent functionality of the wheelchair for the user is the most important task of Veldink4kids.*



*Back side with dynamic back*

## Footrest

Choice of various foot and leg support systems. The continuous footplate and the separate footrests. The footplates can be provided with a heel edge and are height and angle adjustable.



Both continuous footplate and separated footplate are available in fixed and folding version and the separate footplates can be provided with a dynamic spring mounting system.

## Rear wheels

Choice of 22, 24 or 26" Quick Release rear wheels with tyres of your choice. Solid or air tyres in standard or outdoor styles.



**Supervisor foot operation**  
This system can be well used as a parking brake and keeps the hands free.


**Facilitator manual**  
This system can also be used as a parking brake by using the locking system on the handles. (not shown in the figure)

## Brakes

Standard version with for the user Veltrom drum brakes whereby the brake lever can be extended. Beside that, choice from the following supervisor brakes.

- Supervisor foot control.
- Supervisor manual control and manual user control.
- Supervisor foot control and user manual control

## Seats

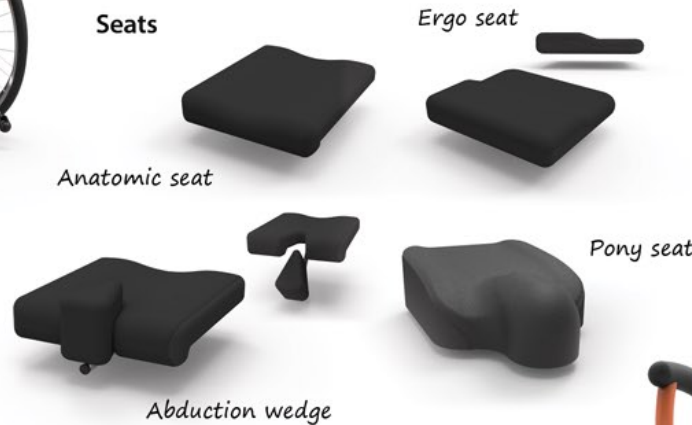

## Backrests




## Headrests

Our different types of headrests ensure optimal position of the head in combination with the right support. The headrest is adjustable in angle, height and depth. Of course we fabricate, if desired, fully customized headrests. By choosing for a stainless steel version, the strength of the position is guaranteed.



## Work trays

The work trays are provided with a double bar to be mounted underneath both armrests. This system is easy to adjust in depth and lock. It is also possible to choose for a swing-away work table, so that the surface can be lifted away.

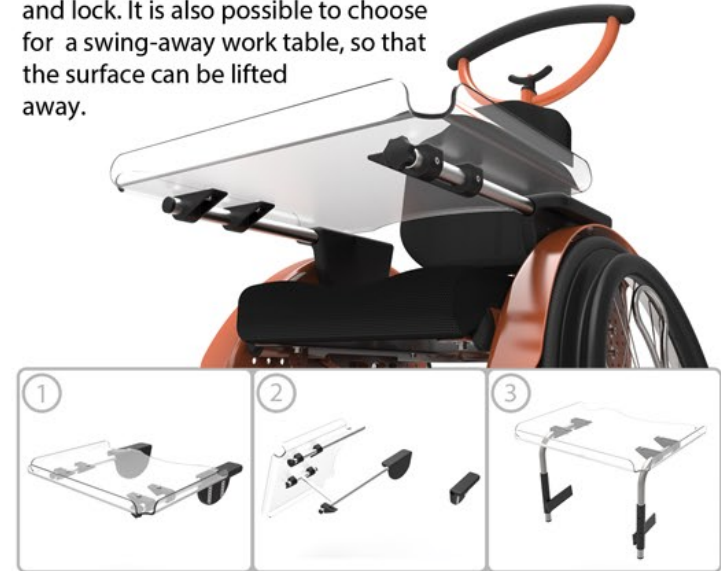

## More options

There are more options possible than displayed in this flyer. Check:

[www.murrays.ie](http://www.murrays.ie)  
[www.murraysmedical.co.uk](http://www.murraysmedical.co.uk)  
[www.veldink4kids.com](http://www.veldink4kids.com)

and the comprehensive option brochures for all the possible options for the Kiddo Classic.