













Specifications

For the Kiddo Space the following data are applicable.

	Seat width 16 - 32 cm		Overall length from 65 cm
	Seat depth 16 - 38 cm		Overall width from 42 cm
	Seat height 34 - 50 cm		Back angle 10 degrees
	Leg length 15 - 40 cm		Camber 7 degrees
	Back height 16 - 40 cm		Kg. Maximum load 50 kg.
	Kg. Overall weight from 7,4 kg included complete seat cushion and back cushion		Rear wheel size 20" 22" 24"
	Kg. Weight without rear wheels from 4,8 kg included complete seat cushion and back cushion		Front wheel size 4" 5"

Standardization

Maximum user weight is 50kg. The Kiddo Space has been tested and approved by TÜV Rheinland for NEN-EN 12182 and NEN-EN 12183 and crash test passed for occupied transportation: ISO7176 -19. Veldink4kids is affiliated with Nationaal Keurmerk Hulpmiddelen (NL)



MURRAY'S

Space
Kiddo by Veldink

The Kiddo Space is a compact wheelchair designed for active users. It is light to self propel and push whilst retaining Veldink Quality and strength. The Kiddo Space is made to measure for a perfect fit with additional growth adjustment available.

Do you want to contact us?



Block C Airton Park,
Airton Rd, Tallaght,
Dublin 24



Eire Tel: 01 855 5696
UK Tel: 028 91 877951



info@murrays.ie
sales@murraysmedical.co.uk

Do you want to follow us?



www.facebook.com/murrays.ie
www.facebook.com/murraysuk



www.twitter.com/murraysmed
www.twitter.com/murraysmeduk



www.murrays.ie
www.murraysmedical.co.uk



For early childhood rehab, framesize from 16cm width



Lightweight wheelchair (from 7kg.)



Growth adjustable in 3 dimensions



Perfectly adjustable to each user's therapeutic and physical demands

Veldink4kids
no limits

Kiddo Space

The Kiddo Space is the ideal children's wheelchair for early childhood, it is manufactured and adapted exactly to the child. The wheelchair grows with the child and is adjustable in width, depth and height.

An innovative children's wheelchair with a contemporary design that is highly manoeuvrable thanks to its compact design and is especially light driving.

Foldable push bar



Fully personalize the wheelchair by the color, wheel protectors and back image designs.

Continuous footplate



The one piece base plate is in height and angle adjustable. The size is determined on the basis of the frame size. The foot plate can be carried out in a fixed and a foldable version. There is also a heel edge mounted version for the footplate to prevent backward movement.

Rear wheels

Choice of 20, 22 or 24" Quick Release rear wheels with tyres of your choice. Solid or air tyres in standard or outdoor styles.

Brakes

Brakes integrated into the mudguards for easy access for the child with optional lock or attendant controlled.

Rear wheel with alu hoop



Outdoor tyre




Seats

We have a wide range of seating options to suit children's varying needs. Each type of seat or backrest cushion is individually customised and made to measure for the child with a wide range of options and additional support accessories.


Abduction wedge



Anatomic seat cushion

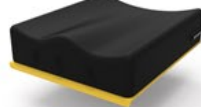

Ergo seat cushion

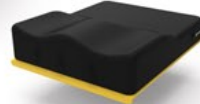

Pony seat



Deep contour back



One piece back



Side support



Spoke protectors

Fully personalise your wheelchair with bespoke spokeguards and trays with virtually unlimited possibilities. We can even incorporate your own images or photos.



Choose your own spoke protectors! www.murrays.ie or www.murraysmedical.co.uk

Headrests

Our different types of headrests ensure optimal position of the head in combination with the right support. The headrest is adjustable in angle, height and depth. Of course we fabricate, if desired, fully customized headrests. By choosing for a stainless steel version, the strength of the position is guaranteed.





Work tray

The work trays are provided with a double bar to be mounted on the fenders. This system is easy to adjust in height and lock on the right height.



Removable work tray, depth and height adjustable. Equipped with double insertion tube for placement in the tubes on the front.

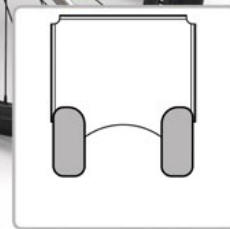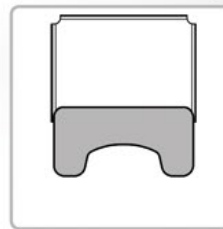

More options

There are more options possible than displayed in this flyer. Check:

www.murrays.ie
www.murraysmedical.co.uk
www.veldink4kids.com

and the comprehensive option brochures for all the possible options for the Kiddo Space.